



POSTOPERATIVE SHOE

Forefoot Off-loading

CHAUSSURE TYPE 1


Original BAROUK *Dr L. Barouk*



- Immediate functional recovery without the use of canes.
- Surgical results protection even with pin in metatarsal phalangeal joint
- Stability and Comfort.
- Bilateral surgery possible with the same functional recovery.

- Usage of clinically tested materials.
- Size: 35/36 — 37/38 — 39/40 — 41/42 — 43/44 — 45/46.
- Specification : Right or Left

TECHNICAL SPECIFICATIONS

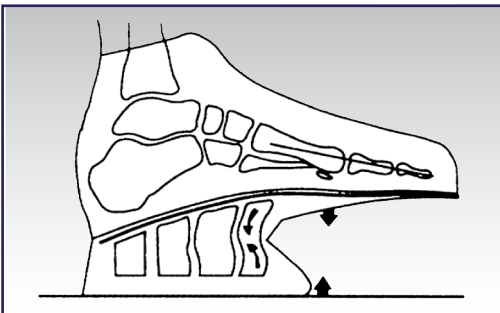

- The talus position of the ankle as well as the dimension of the shoe sole provide tension relief and protection of the fore-foot, even in the event of metatarsophalangeal broaching.
- The heel block is compressible forward allowing a sufficient roll of the foot. The transverse and longitudinal block dimensions ensure the required stability, avoiding backward/forward loss of balance. Working together with the heel block, a shape memory board allows the lack of deformation while using the shoe.

INDICATIONS

- Post operative forefoot surgery (first month).
- Each acute state, inflammatory, painful, oedematous, with the need of a forefoot discharge.
- Continuous pain.

Romans Industrie